


Options

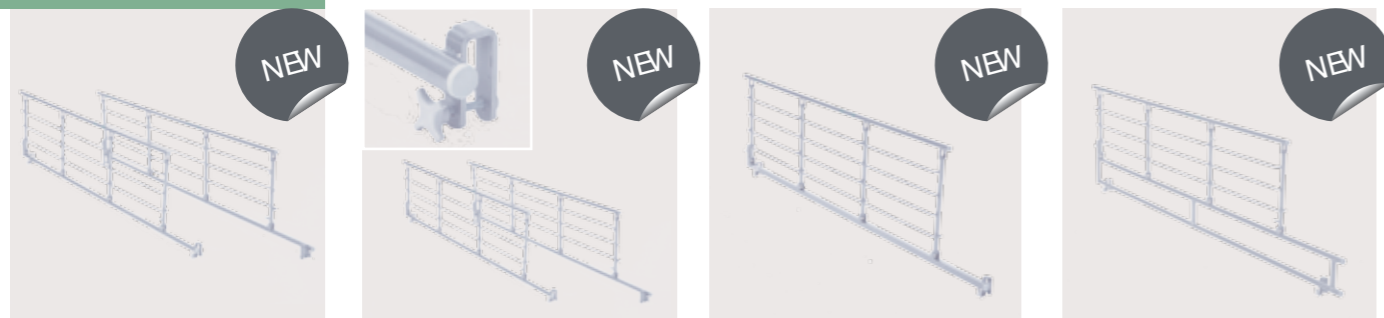

Invacare Scala Basic 2

Foldable steel side rail. Compliant with the new standard for beds IEC 60601-2-52.
Part: 1528930-0154

Invacare Scala Basic Plus 2

Foldable steel side rail with plastic inserts at the mounting brackets. Compliant with the new standard for beds IEC 60601-2-52.
Part: On request

Invacare Scala Medium 2

Foldable steel side rail for deeper mattresses. Compliant with the new standard for beds IEC 60601-2-52.
Part: On request

Invacare Scala Decubi 2

Foldable steel side rail for deep mattresses such as Dynamic systems. Compliant with the new standard for beds IEC 60601-2-52.
Part: On request

Technical data


Etude Plus Width outside: 1020 mm Width inside: 900 mm	Length outside: 2190 mm Length inside: 2070 mm Mattress: 2000 mm	Height adjustable: 330-730 or 400-800 mm	Mattress support dimensions: 850-240-350-560 mm	0 - 70° Backrest angle	0 - 18° Thigh angle

Etude Plus 0 - 25° Legrest lift	11° Tilt option**	75,4 kg	SWL 180 kg / 28.5 stone Max. user weight 145 kg / 23 stone	Heaviest part of the product: Head section 24.6kg

* A reduction in shear and friction forces may contribute to the reduction of pressure ulcers.
** Please note that the tilt function is only suitable for users under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work

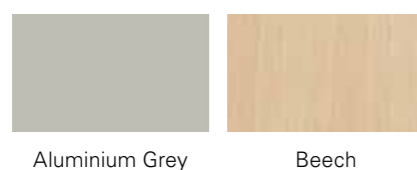

Aluminium Grey Beech

Products, Options and Features can be subject to change.

Warranty

3 years frame
2 years actuators

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk

Invacare Australia Pty Ltd
1 Lenton Place (PO Box 5002)
North Rocks NSW 2151, Australia
Freephone: 1800 460 460
Fax: 02 8839 5353
Email: sales@invacare.com.au
Web: www.invacare.com.au

Invacare New Zealand
4 Westfield Place (PO Box 62-124)
Mt Wellington Auckland, New Zealand
Freephone: 0508 468 222
Freefax: 0508 807 788
Email: sales@invacare.co.nz
Web: www.invacare.co.nz

© 2011 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without prior consultation.

PB5402/1111



NEW

Shear Management Technology

3 year warranty for frame

2 year warranty for actuators

Etude Plus - Leading the way with comfort, ease of use and durability for the community

The **Etude Plus** care bed, with its attractive aesthetics, flexible design and easy handling is ideal for community care.

All the hallmarks and proven features of the popular **Etude**, such as sturdiness, ease of handling and durability, are very much present in this new community bed.

The **Etude Plus** is available with the innovative 'Shear Management Technology' which helps to reduce shear and friction forces. The backrest auto-regression mechanism swings up and back as the bed profiles, reducing the propensity for the client to slide down the bed. This reduction in shear and friction forces may help in the prevention of skin damage and associated pressure ulcer formation. For more information, please see our 'Shear Management Technology' video at www.invacare.co.uk.





Modern and Flexible

The streamlined, discreet design perfectly complements any interior. The bed can be configured with steel or wood side rails with the option of having an extra wood plate to cover the bed end actuator if required. Small details make a big difference to personal comfort: the **Etude Plus** is equipped with auto-regression for increased comfort and reduction in shear and friction*. A lockable hand control is also available for increased safety. The standard, built-in double-brackets are designed to suit any client's transfer needs by offering different positions (entrance height of 400 mm or a lower alternative of 330 mm.)

Easy handling & storage

Portability and storage are key benefits of the **Etude Plus**. Its simplistic design allows easy assembly and dismantling without tools, and the process can be safely and comfortably managed by one person (each section weighs less than 25 kg). The bed ends can be easily attached to the transport unit making the **Etude Plus** lightweight and manoeuvrable for handling and storage.

Highest possible safety levels

The **Etude Plus** conforms mechanically to IEC60601-2-52, the new standard for beds. This equates to the highest degree of safety for both client and carer. The **Etude Plus** is compatible with three different ranges of side rail which all

meet the demands for increased strength required in the new safety standards for beds. For more information on what the new IEC standards mean for you, call Customer Services on 01656 776222.

Options



Etude plus - At home in the home

An emphasis on design ensures the **Etude Plus** blends into home and residential environments perfectly. To further enhance the stylish design features, a wooden panel can be added to cover the bed end actuators.



Turnable Support Handle

The turnable support handle provides good support when moving into a standing position. Part: 1529665-0154



Britt V

Full-length side rails in beech. Compliant with IEC 60601-2-52. Part: 1538458-0101



Line/Line Extended

Full-length aluminium side rails. Compliant with IEC 60601-2-52. Part: 1522786-0125



Height Extender

Height of full-length side rails can be increased by 150 mm. Part: On request

Features and Options



Built-In low version with Double Brackets

Easy and quick attachment of mattress support to bedends in either 400-800 mm or 330-730 mm adjustment height.



Easy Maintenance

IPX6 motors can be dismantled without tools. The cabling between the control unit and the motor is easy to access and replace.



Shear Management Technology

Improved comfort when seated, reduced shear/friction.



Hand Control

Slim, ergonomic hand control. Choose with or without lockable functions (lockable handset available as an option).



Lifting Pole

Light and easily adjustable. Part: 1542093-0154



Transport Kit

For easy storage and enhanced portability. Part: 1543779-0154



Mattress Support Extension

For the ultimate in customisation and personal comfort, an extension of 200 mm is available. Part: 1423994-0154



Cover for wooden side rail

Padded cover for wooden side rails. Cover remains in place when side rail is lowered. Available in standard, padded and net. Part: On request