



Hoist Options



Standard Yoke

- 2 sling attachment points at extremity of curved yoke bar
- Use with general purpose (loop attachment) slings
- Non-absorbent polyurethane padding

Available on: Alto, Sonata & Sonata ELS

Indications: Suitable for the majority of patients

Useful for: All day-to-day patient lifting

Standard Yoke



4 Point Yoke

- 4 sling attachment points spread over a broad area
- Use with general purpose (loop attachment) slings
- Padded vinyl covers with velcro closure

Available on: Alto, Sonata & Sonata ELS

Indications: Bariatric clients, clients requiring a semi-supine lift to maintain upper body control

Useful for: Supine lifting devices eg. Jordan frame

4 Point Yoke



**Innovation in
Patient Handling**

Allegro
CONCEPTS PTY LTD

Hoist Options

The Allegro Concepts range of standard patient lifting hoists includes:

- Alto Lift (SWL 200kg)
- Sonata (SWL 150kg)
- Sonata ELS (SWL 150kg)

These three heavy duty hoists are available with a choice of 6 different sling attachment options.

Achieve quality patient transfers by selecting the fitting which best suits your client or healthcare facility.

Factors which may determine choice of sling attachment include:

- Physical attributes of patients eg. weight, disability
- Level of staff training eg. competence with complex lifting procedures
- Necessity to interface with existing equipment eg. beds, wheelchairs

Tri-Hook Yoke

- 6 Sling attachment points at extremity of curved yoke bar
- Use with general purpose (loop attachment) slings
- Non-absorbent polyurethane padding

Available on: Alto, Sonata & Sonata ELS

Indications: Where the sling prescribed has more than the usual 4 attachment loops

Useful for: Facilities which have historically trained with this configuration

Hi-Lift Yoke

- 2 sling attachment points at extremity of slightly curved yoke bar
- Use with general purpose (loop attachment) slings
- Padded vinyl cover with velcro closure

Available on: Alto, Sonata & Sonata ELS

Indications: Where lifting height of specified hoist is insufficient (increases lifting height by 80mm)

Useful for: Reducing the loss of lifting height when installing an Allegro Universal Weighing Device

Pivot Frame

- 4 Sling attachment points (2 at upper extremities & 2 in lower centre of lifting cradle)
- Use with pivot slings (clip or Eclip attachment)
- Non-absorbent polyurethane padding

Available on: Alto, Sonata & Sonata ELS

Indications: Where patient needs to be tilted for upper body stability, where patient needs to be sat upright to facilitate good positioning

Useful for: Large hospitals with many staff – limited variation in sling attachment options

Powered Pivot Frame

- 4 Sling attachment points (2 at upper extremities & 2 in lower centre of lifting cradle)
- Use with pivot slings (clip or Eclip attachment)
- Non-absorbent polyurethane padding

Available on: Alto & Sonata ELS

Indications: As per regular pivot frame but where manual tilting of the pivot frame is identified as an OH&S risk

Useful for: Preventing inadvertent tilting of pivot frame


Tri-Hook Yoke




Hi-Lift Yoke



Pivot Frame

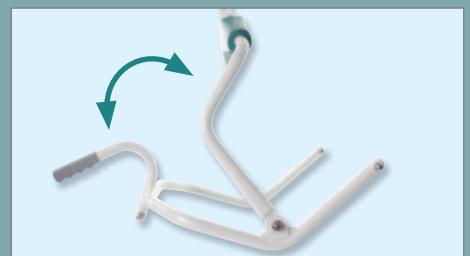

Powered Pivot Frame

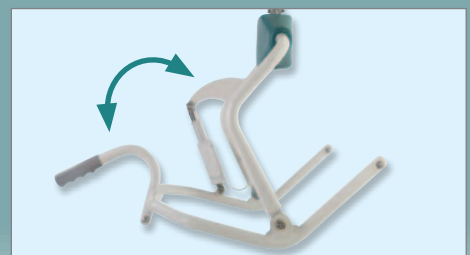

Allegro Concepts Pty Ltd

30-32 James Street Lidcombe NSW 2141

Phone: (02) 9749 7812 Fax (02) 9749 2144

Email: sales@allegroconcepts.com.au