

SALSA 200

Compact Active Standing Aid



Featuring Linak Jumbo Care advanced diagnostic control system
Powered by LA44 European Standards compatible linear actuator

**Innovation in
Patient Handling**

Allegro
CONCEPTS PTY LTD

SALSA 200

SALSA 200: Designed as an active standing aid allowing a single carer to raise a patient into a comfortable secure standing position.

The Salsa 200 is ideal for transferring from bed to chair, or chair to commode. Allegro Concepts standing slings are manufactured to suit the Salsa, allowing superior access for dressing and personal hygiene.

With a small footprint and low overall weight, the Salsa 200 is highly maneuverable. A high capacity safe working load of 200kg means the Salsa is equally at home in an institutional or home environment.



Multiple Grip Positions




Adjustable Kneepad



Removable Footplate



Tempo Handle Kit

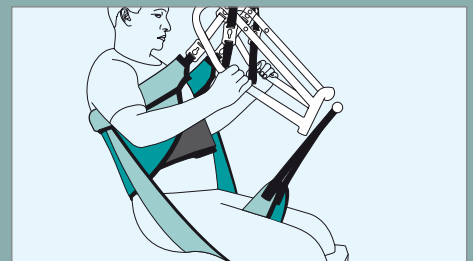

Custom Sitting Sling



Knob Kit for Custom Sitting Sling

Salsa 200 Specifications

Safe Working Load	200kg
Base Length	1010mm
External Base Width Closed	665mm
External Base Width Open	1030mm
Internal Base Width Closed	565mm
Internal Base Width Open	930mm
Base Width Adjustment	Electric
Under Bed Clearance	130mm
Minimum Height Lift Arm	880mm
Maximum Height Lift Arm	1685mm
Under Base Height	55mm
Lifter Weight	51kg
Minimum Storage Height	1190mm
Footplate	Removable
Kneepad	Height Adjustable

The Allegro Concepts range of patient lifting equipment including lifters, slings and accessories has been designed and tested to comply with: AS ISO 10535-2002 Hoists for the transfer of disabled persons – Requirements and test methods.

Allegro Concepts Pty Ltd

30-32 James Street Lidcombe NSW 2141
Phone: (02) 9749 7812 Fax (02) 9749 2144
Email: sales@allegroconcepts.com.au

Allegro
CONCEPTS PTY LTD
www.allegroconcepts.com.au